

iSperm Nexus DX1

English version

Human edition



PLEASE READ THIS MANUAL CAREFULLY BEFORE OPERATION

iSperm

iSperm Nexus DX1 User's Guide - English Version – HUMAN edition

Version: 1

Language: English

Release date: July 2026

© 2026 iSperm Medical



iSperm Medical
www.isperm.com - market@isperm.com

Table of contents

1	Introduction	16
1.1	Product information.....	16
1.2	Product structure and components	16
1.3	Product applicability	17
1.4	Equipment technique parameter	17
2	Equipment description	18
2.1	Front view	18
2.2	Back view	20
2.3	Other components	21
3	Equipment installation	22
3.1	Equipment opening and sealing.....	22
3.2	Select proper position	25
3.3	Power source requirement	25
4	First use	26
4.1	Main screen.....	26
4.2	Settings	27
4.2.1	General.....	27
4.2.2	Analyze.....	32
4.2.3	Reports.....	40
4.2.4	Users	47

iSperm

Intended use



The device is intended for the automated acquisition and processing of microscopic images of in vitro human semen samples.

The device generates structured, image-derived data representing observable characteristics of the samples, to support laboratory-related activities such as method standardization, workflow consistency control, and internal reference purposes.

The device does not provide medical interpretation, is not intended to establish or exclude any clinical condition, and shall not be used as an independent basis for diagnostic or therapeutic decision-making.

The outputs generated by the device are non-decisional and auxiliary in nature, and are intended to be used by qualified professionals in conjunction with other laboratory information and relevant contextual data.

How to use this manual

It's greatly appreciated you become the user of Nexus Dx1.

To get the best effect, you must be familiar with our equipment and their functions before performing the clinical diagnostic tests.

This manual is the operating guidance of Nexus Dx1 semen analyser, which includes the contents of equipment installation, daily test, quality control and maintenance.

Please, reserve all the packing materials for the use of deposit, transportation and maintenance in the future.

You may contact your distributor if there is any problem.

Warning:



Indications must be obeyed, to ensure the normal operation of equipment and correct and real test results.

Symbols on the equipment



This means that the labelled item could lead to personal injury and/or damage to the analyser. The symbol is labelled beside the power outlet and some external interface.



The symbols for "SERIAL NUMBER", The serial number shall be after or below the symbol, adjacent to it.



The symbol indicates the manufacturer and its address, after which are shown its name and address.



The symbol indicates the date of manufacture, after which are shown its manufacture date.



The symbol means that the environment of equipment must be damp proof during transport, and equipment must be kept in a dry environment.



This means that equipment should handle with care during transportation, so as not to damage it.



The symbol means the equipment packaged should not be upended at any time.



This symbol indicates that when the end-user wishes to discard this product, it must be sent to separate collection facilities for recovery and recycling. By separating this product from other household-type waste, the volume of waste sent to incinerators or land-fills will be reduced and natural resources will thus be conserved.

Warning and safety indications

This equipment is only used for *in vitro* diagnosis, please carefully read the following warnings before the usage, and strictly obeys them.



Please carefully read the following announcements before using this equipment.

Cut off the power source at once if there are abnormal smells or smokes and pull out the power plug from the outlet. At that moment, you must immediately propose test applications to the distributors and the agents of our company. If it is continuously used under this condition, it may cause fire, electric shock or personnel casualty.

Neither water and reagent nor metal pieces of staple and pin should be dropped into the equipment; otherwise it would lead to short circuit or smoke and fire. If there is something abnormal, cut off the power source at once, and pull out the power plug from the outlet. At that moment, you must immediately propose test applications to the distributors and the agents of our company.

Operators don't touch the electric circuit inside the equipment, especially wet hands, which may be easier to get electric shock.

Rubber gloves must be put on while maintaining and checking the equipment, and prescriptive tools and spare parts should be used. After finishing the operation, please use disinfecting liquid soap to wash hands.

The used tubes and other equipment should be properly disposed as the medical wastes or infective wastes. If the person is stained with blood, there is possibility of being infected by pathogens.

Power voltage, connection and earthing

Be sure not to insert the plug into the electrical outlet of alternating current over 220v. Otherwise, it may cause fire or electrical shock.

When the equipment is installed, three-core power cable accompanied with the machine must be used, and ensure it is well earthing. Otherwise, it may cause fire or electrical shock.

Be sure not to damage the isolation protection peel of power line. Do not pull the power line with exertion or hang heavy things on it. Otherwise, it may cause broken circuit or short circuit, thus lead to electrical shock or fire.

Power source must be cut off when connect the peripheral equipment. Otherwise, it may cause electrical shock or troubles.

Shelf-Life

Once the device is put into operation, there is a shelf life of 6 years. The 6-year device life-time limit is applicable if the final user chooses NOT to install software updates that may be developed by the constant evolution of technology and incorporation of new clinical tests.

1 Introduction

1.1 Product information

Product name: Nexus Dx1

Model: v.1.0

Product property:

Nexus Dx1 is an automate system intended for the automated acquisition and processing of microscopic images of in vitro human semen samples.:

- Concentration: Calculates the total sperm number of a semen sample and the sperm concentration (M/mL).
- Motility: Calculates the total percentage of motile sperm.
- Morphology: Calculates the total percentage of morphology pattern sperm.
- Round cells pattern: Calculates the concentration of round cells like Non-sperm particles.

After the analysis, the system sends all the results, videos, images and sessions to the Laboratory Information System (LIS), where is possible to check the analysis or print reports.

1.2 Product structure and components

Nexus Dx1 is internally composed by:

Optical system

Motorized stage

Warm stage

Digital camera

Computer control system and software

1.3 Product applicability

Nexus Dx1 is intended for the automated acquisition and processing of microscopic images of in vitro human semen samples.

1.4 Equipment technique parameter

Light source: LED light

Temperature control: 35 - 40 °C (37°C as default)

Interface: RJ45 protocol

Display: Internal displayer 7 inches

Print: External printer

Work environment: 0°C to 40°C; relative humidity<85%

Storage environment: -10°C to 40°C; relative humidity<85%

Weight: 7kg

Overall dimensions: 283mm(L)207mm(W)x244mm(H)

Power supply: A.C. 110V/240V, 50Hz/60Hz

Power input: 12V 5A

2 Equipment description


2.1 Front view




Figure 1. Nexus Dx1 front view

1. **Indicating lamp (left):** Power and function lamp.


- Always on if power switch is on.

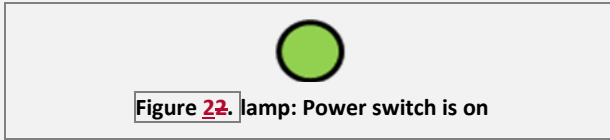


- Off if power switch on the rear panel is off.

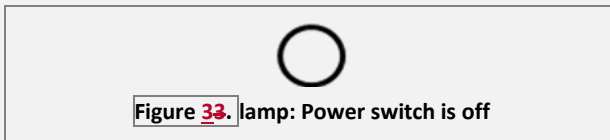


- Green: Always on if the system is ok.

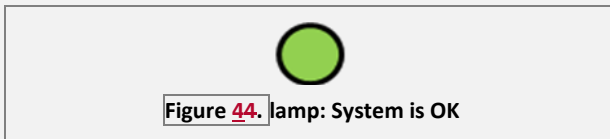

- Red: Be careful. Warning or problem inside the hardware system.

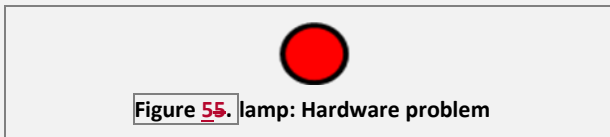

- Blue: Be careful. Warning or problem inside the software system.

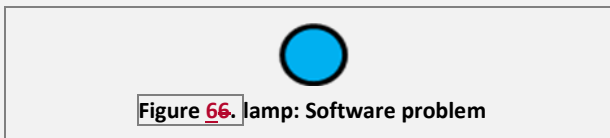

2. **Sleep/Wake-up Button (Right)**

Press this button to put the screen into sleep mode. When the screen enters sleep state manually or automatically, press this button to wake the screen up.

2.2 Back view

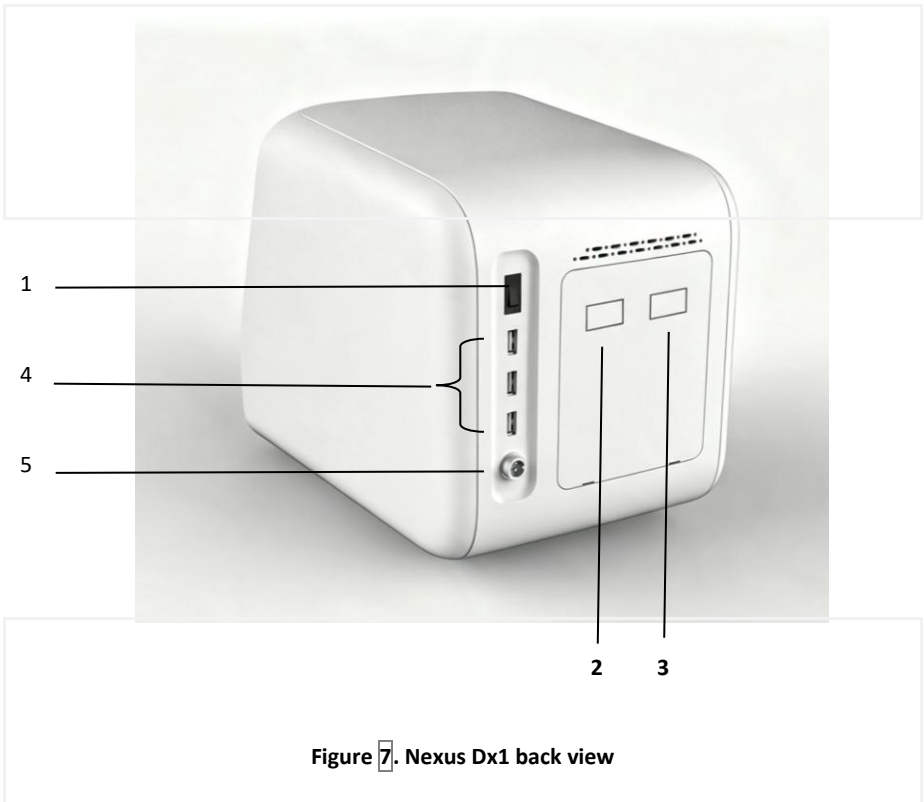


Figure 7. Nexus Dx1 back view

1. Power switch
2. PC Monitor connection (Screen Port)
3. PC software connection (Data port)
4. 3 USB ports. Connect keyboard, mouse, barcode scanner, USB disk
5. Power cable

2.3 Other components

Nexus Dx1 also includes the following components:

1. Barcode scanner wired
2. Keyboard + mouse
3. USB flash drive (for APP update, Export Reports), FAT32 Only, Device-Specific
4. Screen port cable (Mirco USB cable)
5. Data port cable (Network RJ45 cable)
6. Power Adapter (Input: 100-240V 50/60Hz 1.4A, Output: +12V 5A 60W)

This system is also compatible with any USB keyboard and mouse which is compatible with PC.

3 Equipment installation

3.1 Equipment opening and sealing

1. Open the transportation box.

Password Lock User Manual



① Initial password 000, push the square button to  to open the password lock and pull out the red seal pin.

② Password change method:
Push the square button to  and hold, rotate the password lock numbers, adjust the numbers, then release the square button.

Note: (Please remember the password after setting it)

Figure 8. Transportation box lock

2. Take off the packing materials and take the equipment out of the package. Be careful: Nexus Dx1 is a heavy equipment. Preserve the packing box and materials well provided that you may need to repack the equipment in the future.



Figure 9. Content of the packing box

3. Check the articles inside the package box, and confirm the following stuffs should be included:
 - Nexus Dx1 Analyzer
 - User's manual
 - Packing list
 - Warranty certificate
 - Power line and AC adapter
 - Keyboard and mouse
 - Barcode scanner
 - RJ45 adapter USB to Ethernet cable
 - Data port cable
 - Screen port cable
 - USB flash drive (FAT32 formatted)

Note: The product accessories should be in accordance with the packing list. Please contact with the distributor when the accessories are not correct.

3.2 Select proper position

Find a place without direct sunshine in your working place. The surface of the selected working table should be flat, with enough space for placing SCA SCOPE and monitor. The anterior edge of equipment should be near to the edge of working table. Avoid comparatively large shake of working table.

Note: The work environment of equipment is that temperature is between 0°C-40°C and relative humidity is below 85%.

To guarantee the normal working of equipment, the following places are forbidden to be put on:

- Places with extreme changes of temperature.
- Places near a centrifuge or shaking machines.
- Particularly hot or cold places.
- Places full of dust.
- Places nearer to the electromagnetic devices that produce magnetic field.
- Places with direct sunshine.

3.3 Power source requirement

A.C.110V/240V

50Hz/60Hz

60VA

4 First use

4.1 Main screen

1. Switch on the Nexus Dx1.

After starting the system, the main screen (7" Built-in LCD Touch Screen) will appear:

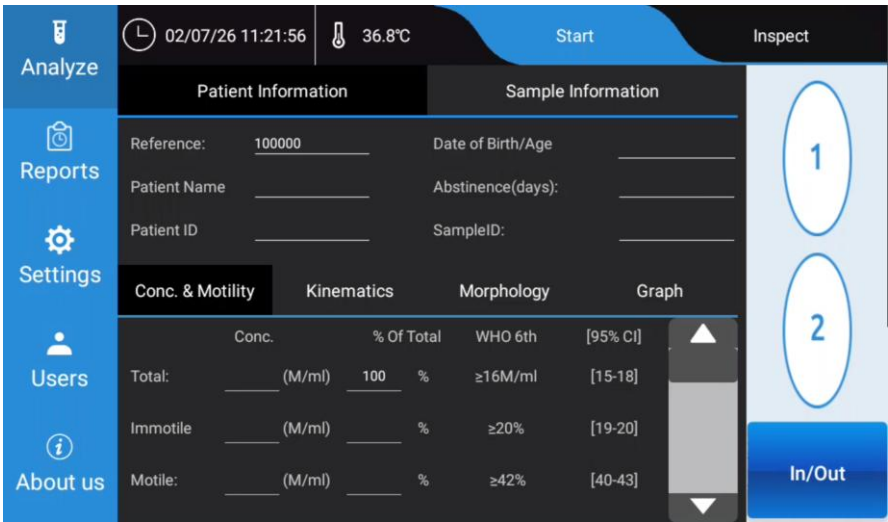


Figure 10. Main screen

Options available:



Language selector: Change the system’s language.



Power options:

- Shut down: Turn off the system.
- Restart: Turn off the system and restart it automatically.



Settings: Access to the system’s configuration options.

4.2 Settings

Click the “Settings” button to enter to the system’s configuration menu. The “Settings” button is available at the main screen as well as the “Analysis” screen.

4.2.1 General

In this section, there is the option of Microscopic Fields, setting to adjust the number of imaging fields per chamber. Available selections are 1, 2, or 3. When you pick a value, the system will analyze that many fields within each sample chamber, then calculate an average value from all captured fields.

Using more fields delivers higher test accuracy, though it will take longer to complete analysis.

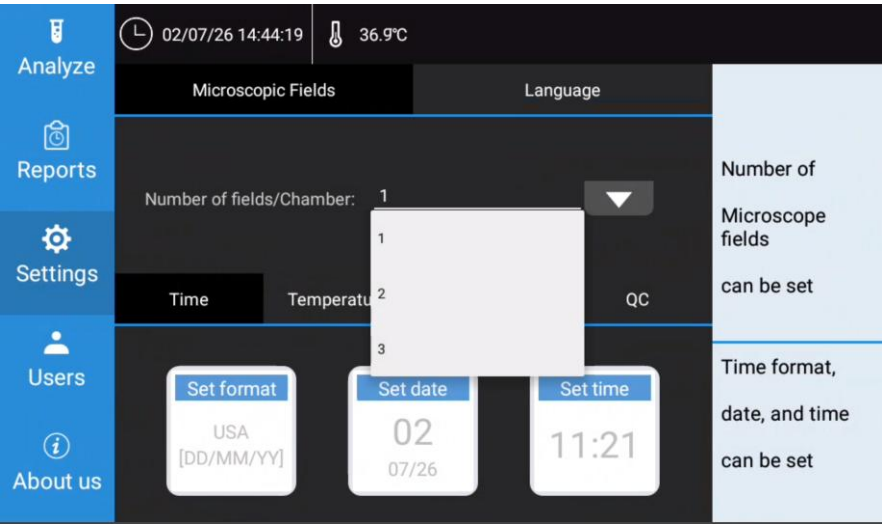


Figure Figure 1032. Settings menu/ Language
1144. Settings menu/ Microscopic Fields

In this section, there is the option of Language selector: Change the system’s language.

A total of twelve different languages are available for selection in the dropdown menu.

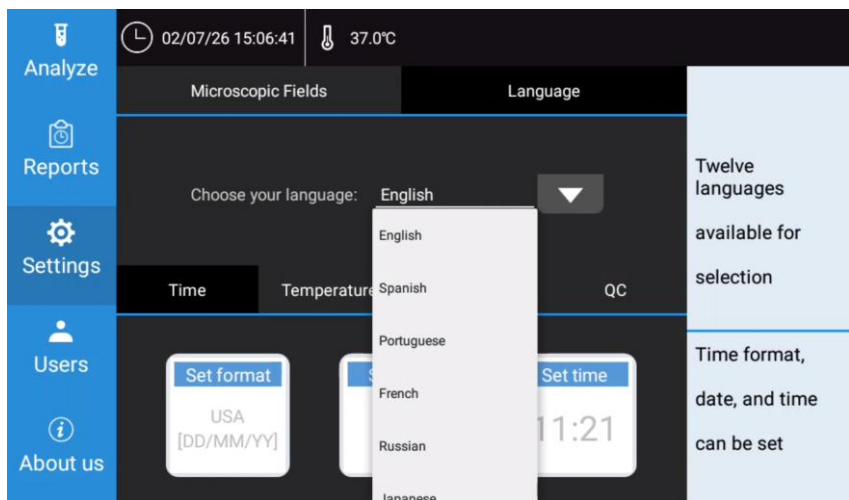


Figure 1212. Settings menu/ Language

In this section, there is the option of Time. Time format, date, and time can be set.


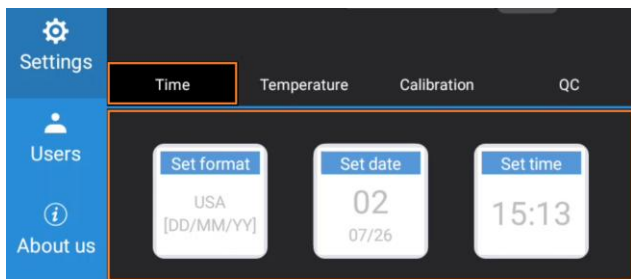


Figure 1313. Settings menu/ Time

In time format section, there are two formats available: USA [MM/DD/YY] and Europe [DD/MM/YY].

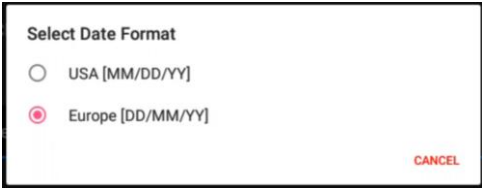

Select Date Format

☐ USA [MM/DD/YY]

☒ Europe [DD/MM/YY]

CANCEL

Figure 1444. Settings menu/ Time / format

In date section, you can set the system’s current date.



Date Selection

Jun	01	2025
Jul	02	2026
Aug	03	2027

CANCEL CONFIRM

Figure 1545. Settings menu/ Time / date

In time section, you can set the system’s current time.

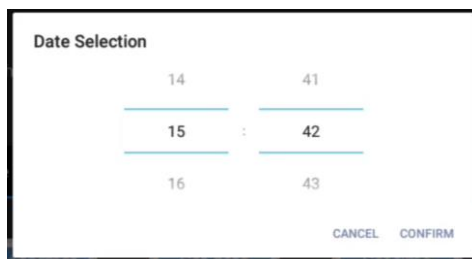


Figure 1616. Settings menu/ Time / time

In Temperature section, the temperature of the platform can be set, the range is 35°C ~ 40°C.

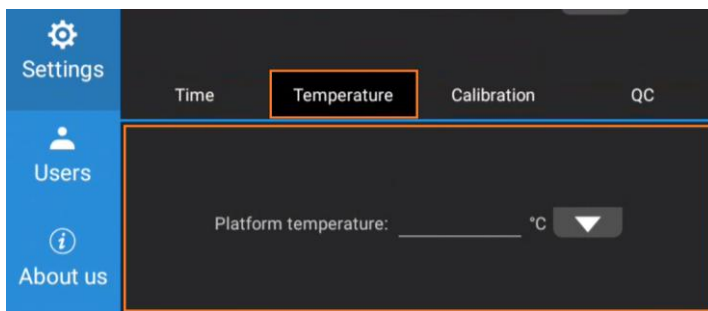


Figure 1717. Settings menu/ Temperature

In Calibration section, the calibration factor of concentration and motility can be set. More information at Calibration guide.

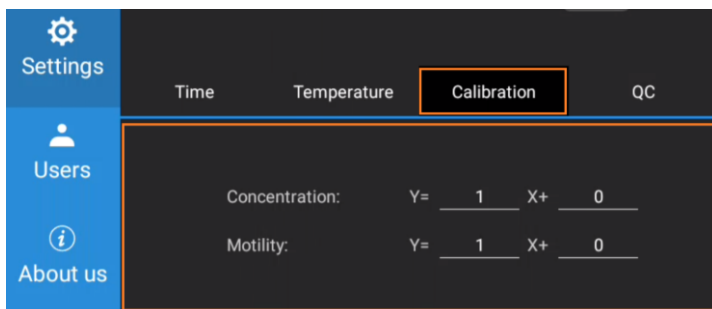


Figure ~~1818~~. Settings menu/ Calibration

In QC section, the QC test using Latex beads can be performed. More information at Calibration guide.

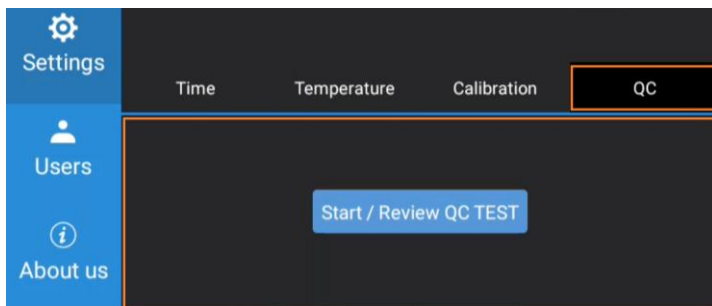


Figure ~~1919~~. Settings menu/ QC

4.2.2 Analyze

4.2.2.1 Overview

The Analyze function is used to inspect a semen sample image, select the chamber to be tested, enter optional patient and sample information, and start an automated sperm quality analysis.

Nexus DX1 is designed for a simplified workflow. The operator may enter patient and sample information before testing, or leave the information fields blank. If no information is entered, the system automatically assigns a unique **Reference** number to the test record. The analysis result can later be found in the **Reports** section by using this reference number.

The analysis result is displayed in four result tabs at the lower part of the screen:

- **Conc. & Motility**
- **Kinematics**
- **Morphology**
- **Graph**

Note: Nexus DX1 provides image-derived analysis data for laboratory use. Final clinical interpretation should be made by qualified professionals according to local laboratory procedures and applicable clinical guidelines.

4.2.2.2 Analyze Screen

From the left navigation bar, tap **Analyze** to enter the analysis function.

The Analyze screen contains the following main areas:

Area	Description
Top status bar	Displays system time and heating stage temperature.
Start / Inspect tabs	Used to inspect the sample and start a test.
Patient Information	Used to enter or confirm patient-related information.
Sample Information	Used to enter or confirm sample-related information.
Chamber selection	Used to select chamber 1 or 2 .
Result tabs	Displays Conc. & Motility, Kinematics, Morphology, and Graph results.

Analyze

02/07/26 11:21:56 36.8°C Start Inspect

Patient Information

Reference: 100000 Date of Birth/Age

Patient Name Abstinence(days):

Patient ID SampleID:

Sample Information

Conc. & Motility Kinematics Morphology Graph

	Conc.	% Of Total	WHO 6th	[95% CI]
Total:	(M/ml)	100 %	≥16M/ml	[15-18]
Immotile	(M/ml)	%	≥20%	[19-20]
Motile:	(M/ml)	%	≥42%	[40-43]

1

2

In/Out

Figure 20. Main screen/ Analyze

4.2.2.3 Inspect the Sample

Before starting a test, the operator should inspect the semen sample image.

1. Tap **Analyze** on the left navigation bar.
2. Tap **Inspect** on the upper-right area of the screen.
3. Select chamber **1** or chamber **2**.
4. Observe the microscopic image of the semen sample.
5. Confirm that the selected chamber contains a suitable sample image for analysis.

Use the chamber that contains the sample image you want to analyze. The test will be performed on the microscopic image from the chamber selected during the Inspect step.

If the sample image is not suitable, prepare or load the sample again according to the laboratory procedure before starting the test.

4.2.2.4 Enter Patient Information

After inspection, tap **Start** to enter or confirm patient and sample information.

In the **Patient Information** section, the operator may enter:

- **Reference**
- **Patient Name**
- **Patient ID**
- **Date of Birth / Age**
- **Abstinence (days)**
- **Sample ID**

The **Reference** field is automatically generated by the system. The operator may use this reference number to locate the result in the **Reports** section after the test.

Patient information is optional. If the fields are left blank, the system can still start the analysis and save the result with the automatically assigned reference number.

The screenshot shows a handheld device interface for patient information entry. At the top, there is a status bar with a clock icon, the date and time '02/07/26 16:17:54', a temperature icon, and the value '36.7°C'. Below this is a navigation bar with 'Analyze' (selected), 'Reports', and a settings icon. The main area is divided into two sections: 'Patient Information' and 'Sample Information'. The 'Patient Information' section contains fields for 'Reference:' (pre-filled with '100000'), 'Patient Name', and 'Patient ID'. The 'Sample Information' section contains fields for 'Date of Birth/Age', 'Abstinence(days):', and 'SampleID:'. A large blue button with the number '1' is visible on the right side of the form. At the bottom, there is a QWERTY keyboard with a 'Next' button.

Figure 21. Patient Information

4.2.2.5 Enter Sample Information

In the **Sample Information** section, the operator may enter:

- **Liquefaction time**
- **Liquefaction treatment**
- **Visual appearance**
- **Viscosity**
- **PH**
- **Ejaculate Volume**
- **Total Sperm Number (system calculate)**

Sample information is optional. However, entering complete information is recommended when the result will be reviewed, printed, exported, or linked to external laboratory records.

The screenshot shows a medical device interface with a dark theme. On the left is a vertical sidebar with three icons and labels: a test tube icon for 'Analyze', a clipboard icon for 'Reports', and a gear icon for 'Settings'. The main display area is divided into two tabs: 'Patient Information' and 'Sample Information', with the latter being active. At the top of the main area, there is a status bar showing a clock icon, the date and time '02/07/26 16:19:25', a temperature icon, '37.0°C', and two buttons: 'Start' and 'Inspect'. Below the tabs, the 'Sample Information' section contains two input fields: 'Liquefaction time' with a text input and '(minutes)' label, and 'Liquefaction treatment' with a dropdown menu showing 'Yes'. To the right of these fields is a vertical stack of three buttons: a blue button with a white '1' inside a circle, a grey button, and a grey button with a downward arrow. At the bottom of the screen is a numeric keypad with buttons for '-', '+', '.', '1', '2', '3', '*', '/', ',', '4', '5', '6', '7', '8', '9', '(', ')', '=', and a 'Next' button. A '0' button and a '#' button are also present at the bottom.

Figure 22. Sample Information

4.2.2.6 Start the Test

After confirming the selected chamber and entering any required information:

1. Tap **Start**.
2. A confirmation message appears:
“Is all information input completed? Start the test?”
3. Tap **YES** to start the analysis.
4. Tap **NO** to return to the information screen and check or edit the information.

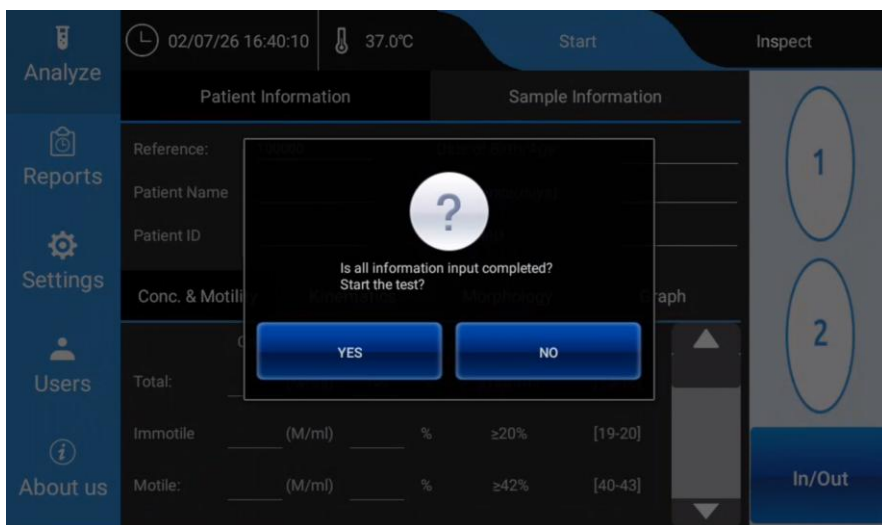


Figure 23. Start test confirmation dialog

When **YES** is selected, Nexus DX1 starts automated analysis of the semen microscopic image in the selected chamber.

Do not remove the sample, switch chambers, or interrupt the instrument during the test unless instructed by the system.

4.2.2.7 View Analysis Results

After the test is completed, the result is shown in the lower part of the Analyze screen. The operator can switch between the four result tabs.

7.1 Conc. & Motility

The **Conc. & Motility** tab displays concentration and motility-related results.

Typical displayed information includes:

- Total (sperm concentration)
- Immotile (sperm)
- Motile (sperm)
- PR (Progressive)
- Rapid-PR (Rapid progressive)
- NP (Non Progressive)
- Round cells (Round cells like Non-sperm particles)

The values shown in this tab are used to review the overall concentration and motility profile of the sample.

7.2 Kinematics

The **Kinematics** tab displays sperm movement parameters calculated from the analyzed image sequence.

Depending on system configuration, kinematic parameters may include:

- DCL (Distance along the curvilinear path)
- DSL (Distance along the straight-line path)
- DAP (Distance along the average path)
- VCL (velocity along the curvilinear path)
- VSL (velocity along the straight-line path)
- VAP (velocity along the average path)
- ALH (the amplitude of the lateral displacement of the head)
- MAD (mean angular displacement (degrees))
- BCF (beat-cross frequency (Hz))
- LIN (linearity of the curvilinear path)
- STR (straightness: linearity of the average path (straight-line velocity / average path velocity))
- WOB (wobble: a measure of oscillation of the curvilinear path about the average path (average path velocity / curvilinear velocity))

These parameters help describe sperm movement pattern and trajectory characteristics.

7.3 Morphology

The **Morphology** tab displays morphology-related analysis results.


Depending on the enabled analysis configuration, the system may provide information related to:

- Normal (sperm)
- Abnormal (sperm)
- Head defects (Head feature distribution)
- Midpiece defects (Midpiece feature distribution)
- Tail defects (Tail feature distribution)
- CR (Cytoplasmic residue)
- TZI (Teratozoospermia index)
- SDI (sperm deformity index)

Morphology results should be reviewed according to the laboratory's validated procedures.

7.4 Graph

The **Graph** tab displays graphical analysis information. This view helps the operator review the sample result visually and include motility chart, morphology chart.

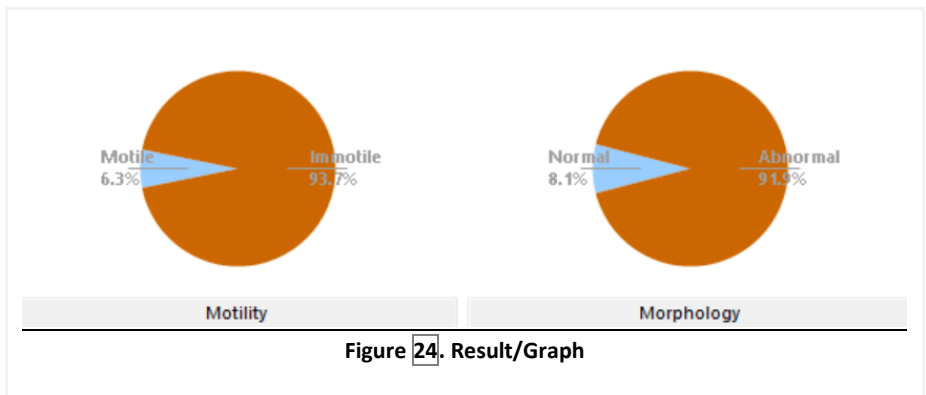

4.2.2.8 Retrieve the Result from Reports

Each completed test is saved by the system. If patient or sample information was not entered, use the automatically generated **Reference** number to identify the record.

To retrieve a result:

1. Tap **Reports** on the left navigation bar.
2. Locate the test record by reference number, patient information, sample information, or test time.
3. Open the record to view, print, export, or review the report according to the available report functions.

4.2.2.9 Recommended Operating Notes

- Inspect the chamber before starting the test.
- Confirm that the correct chamber, **1** or **2**, is selected.
- Enter patient and sample information when traceability is required.
- If information is not entered, record or note the automatically generated reference number.
- Do not disturb the instrument during analysis.
- Review all result tabs before saving, printing, or exporting the report.
- Follow the laboratory's sample preparation and quality control procedures.

4.2.2.10 Troubleshooting

Situation	Possible Cause	Recommended Action
The operator cannot find the result after testing.	Patient or sample information was not entered.	Search in Reports using the automatically generated reference number or test time.
The wrong chamber was analyzed.	Chamber selection was not confirmed during Inspect.	Repeat the analysis with the correct chamber selected.
The confirmation dialog appears before the operator is ready.	Start was tapped before information entry was completed.	Tap NO , complete or confirm the information, then tap Start again.
Result fields are incomplete or unexpected.	Sample image, chamber loading, or analysis conditions may not be	Inspect the sample image again, prepare a new chamber if needed, and repeat the test according to

Situation	Possible Cause	Recommended Action
	suitable.	laboratory procedure.

4.2.3 Reports

Configuration options for review the reports data of Nexus DX1 . This module is always included in Nexus DX1 .

The Reports function is used to review, open, delete, and export completed test results. Each completed analysis is saved as a report record. The operator can review the digital result data, open the full report, play back the analysis video, delete selected records, or export reports to the USB flash disk supplied with the instrument.

The Reports screen includes two main viewing modes:

- Digital Results: Displays the report list and numerical result details.
- Video Playback: Displays the recorded microscopic video of the selected test.

The right-side function buttons include:

- Report List
- Full Report
- Delete
- Export

The screenshot displays the 'Reports' screen of the Nexus DX1 instrument. The interface is divided into a left sidebar, a main content area, and a right sidebar with function buttons.

Left Sidebar (Navigation):

- Analyze (Icon: Test tube)
- Reports (Icon: Document with magnifying glass, currently selected)
- Settings (Icon: Gear)
- Users (Icon: Person)
- About us (Icon: Information 'i')

Main Content Area:

At the top, there is a header bar with the following information:

- Time: 12/02/26 08:41:01
- Temperature: 36.8°C
- Mode: Video Playback (selected) / Digital Results

Below the header is a table of test results:

Ref.	Date	Time	Patient Name	Patient ID	CONC. Total(M/ml)	MOT Motile(%)	MOR Normal(%)
100004	2026-02-12	08:40	Benny	yy123	20.430	68.8	35.910
100003	2026-02-11	16:11	David	kk123	20.430	68.8	35.714
100002	2026-02-09	21:06	Michael	Yy 123	20.430	68.8	26.588

Right Sidebar (Function Buttons):

- Report List
- Full Report
- Delete
- Export

Figure 25. Reports Main Screen

4.2.3.1 Open the Reports Function

1. Tap Reports on the left navigation bar.
2. The system enters the Reports function area.
3. Tap Report List if the report list is not already displayed.
4. Confirm that Digital Results is selected at the top of the screen.

The main area displays all saved test results in a list.

4.2.3.2 Report List

The Report List displays saved test records. Each row represents one completed test.

Typical columns include:

- Ref.
- Date
- Time
- Patient Name
- Patient ID
- CONC. Total (M/ml)
- MOT Motile (%)
- MOR Normal (%)

To select a report, tap one row in the list. The selected row is used for Full Report, Delete, and Video Playback operations.

4.2.3.3 Open a Full Report

To open the full report of a selected result:

1. Enter Reports.
2. Tap Report List.

3. Select one report row from the list.
4. Tap Full Report on the right side.

The screen changes from the report list to the selected report's digital result view.

The screenshot displays the 'Full Report - Digital Results' screen. On the left is a blue sidebar with icons for Analyze, Reports, Settings, Users, and About us. The top header shows a clock icon, the date and time '12/02/26 08:41:07', a temperature icon with '36.8°C', and tabs for 'Video Playback' and 'Digital Results'. Below the header, there are two main sections: 'Patient Information' and 'Sample Information'. The 'Patient Information' section includes fields for Reference (100004), Patient Name (Benny), and Patient ID (yy123). The 'Sample Information' section includes Date of Birth/Age (1995-02-12), Abstinence(days) (6), and SampleID (ff123). Below these is a tabbed interface with four tabs: 'Conc. & Motility' (selected), 'Kinematics', 'Morphology', and 'Graph'. The 'Conc. & Motility' tab displays a table with columns: Conc., % Of Total, WHO 6th, and [95% CI]. The table has three rows: 'Total', 'Immotile(D)', and 'Motile'. To the right of the table is a vertical scrollbar. On the far right, there is a light blue sidebar with four blue buttons: 'Report List', 'Full Report', 'Delete', and 'Export'.

	Conc.	% Of Total	WHO 6th	[95% CI]
Total:	20.430 (M/ml)	100.0 %	≥16M/ml	[15-18]
Immotile(D)	6.374 (M/ml)	31.2 %	≥20%	[19-20]
Motile:	14.056 (M/ml)	68.8 %	≥42%	[40-43]

Figure 26. Full Report - Digital Results

4.2.3.4 Review Digital Result Data

In the full report view, the patient and sample information is shown in the upper part of the screen. The digital result data is shown in the lower part of the screen.

The operator can switch between the following result tabs:

- **Conc. & Motility**
- **Kinematics**
- **Morphology**
- **Graph**

Tap each tab to review the corresponding result category for the selected test.

4.1 Conc. & Motility

The **Conc. & Motility** tab displays concentration and motility results, such as total concentration, immotile sperm, motile sperm, percentage of total, WHO 6th Edition reference information, and confidence interval information where available.

4.2 Kinematics

The **Kinematics** tab displays sperm motion parameters calculated by the system. These values describe the movement characteristics of the analyzed sperm cells.

4.3 Morphology

The **Morphology** tab displays morphology-related analysis results for the selected test.

4.4 Graph

The **Graph** tab displays graphical summaries of the selected result, such as motility and morphology charts.

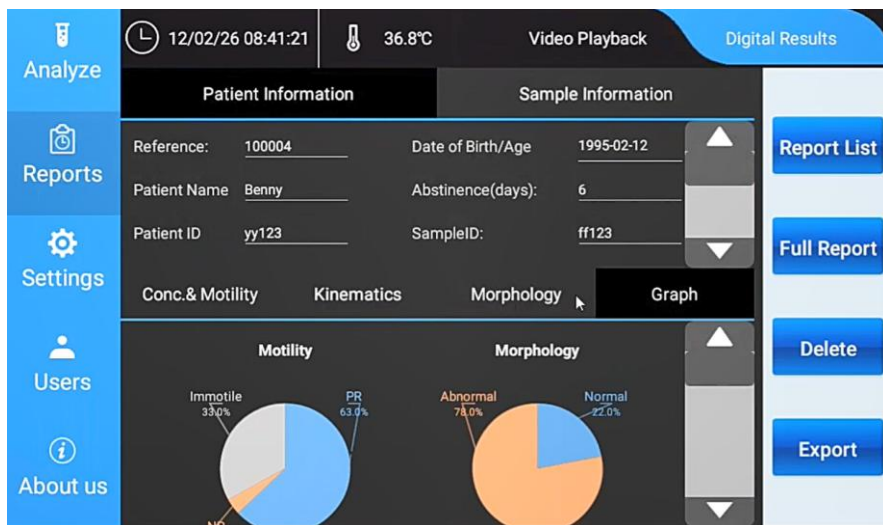


Figure 27. Full Report - Graph tab

4.2.3.5 Video Playback

The operator can review the recorded microscopic video for a selected test.

To open video playback:

1. Open **Reports**.
2. Select a report from the report list, or open the report in **Full Report** view.
3. Tap **Video Playback** at the top of the screen.

The screen displays the microscopic video recorded during analysis.

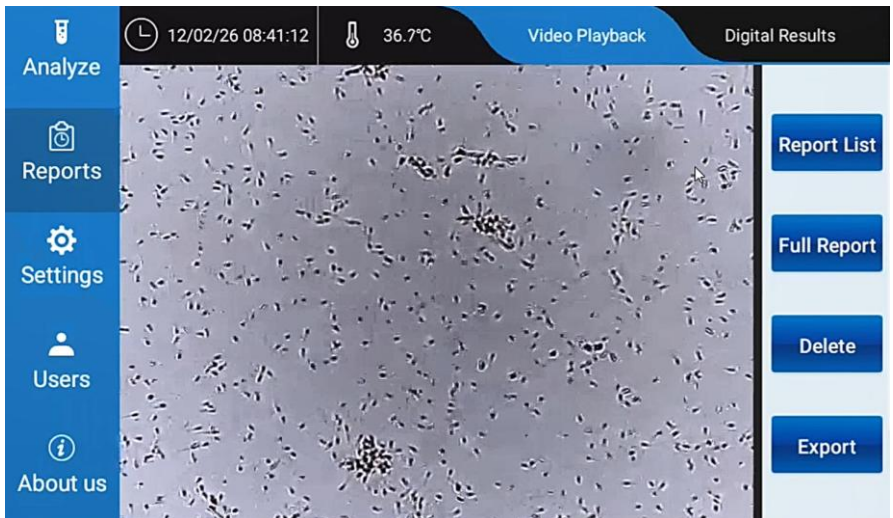


Figure 28. Video Playback

Use video playback to visually review the analyzed sample image. This function is helpful for record review, training, and troubleshooting.

Delete a Report The operator can review the recorded microscopic video for a selected test.

4.2.3.6 Delete a Report

To delete a saved report:

1. Tap **Reports** on the left navigation bar.
2. Tap **Report List**.
3. Select the report row that should be deleted.
4. Tap **Delete** on the right side.
5. Confirm the deletion if a confirmation message appears.

Caution: Deleted reports may not be recoverable. Confirm that the correct report is selected before deleting it.

4.2.3.7 Export Reports to USB Flash Disk

The **Export** function downloads saved report files to the USB flash disk supplied with the instrument.

To export reports:

4. Insert the USB flash disk supplied with the instrument into a USB port on the back of Nexus DX1.
5. Tap **Reports** on the left navigation bar.
6. Tap **Export** on the right side.
7. When the password prompt appears, enter:

000000 (default)
8. Confirm the export operation.
9. Wait until the system finishes downloading the report files to the USB flash disk.
10. Remove the USB flash disk only after the export is completed.

Note: Use the USB flash disk supplied with the instrument when exporting reports. Do not remove the USB flash disk during export.

4.2.3.8 Recommended Operating Notes

- Use **Report List** to locate saved test results.
- Select the correct row before tapping **Full Report**, **Delete**, or **Video Playback**.
- Use **Digital Results** to review numerical result data.
- Use **Video Playback** to review the recorded microscopic video.
- Insert the supplied USB flash disk before using **Export**.
- Enter password **000000 (default)** when prompted during export.
- Confirm that export is completed before removing the USB flash disk.

4.2.3.9 Troubleshooting

Situation	Possible Cause	Recommended Action
The expected report is not visible in the list.	The result may be listed under its automatically assigned reference number or another test time.	Check the Ref. , Date , and Time columns.
Full Report does not open.	No report row is selected.	Select one report row, then tap Full Report again.
The wrong report was opened.	A different row was selected in the report list.	Return to Report List , select the correct row, and tap Full Report .
Video playback is not shown.	The operator is still in Digital Results view.	Tap Video Playback at the top of the screen.
Export does not start.	USB flash disk is not inserted or not recognized.	Insert the supplied USB flash disk into the USB port on the back of the instrument, then try again.
Password is requested during export.	Export is password protected.	Enter 000000 (default) .

4.2.4 Users

The **Users** function is used to manage the instrument password and edit local instrument information. This information is stored and displayed on the instrument only.

In the current software version, the information edited in this section is used for local reference and viewing only. It does not provide additional workflow, reporting, or remote management functions.

4.2.4.1 Open the Users Function

1. Tap **Users** on the left navigation bar.
2. The system enters the **Password Management** screen.
3. The right side of the screen displays brief notes for the current input area.

Analyze

Reports

Settings

Users

About us

12/02/26 08:42:06

36.8°C

Password Management

Old password:

New password:

Confirm

Company

Center

Laboratory

Technician

Company:

Enter 6 to 15
Digital numbers
as password

Enter or select
Company Title

Figure 29. Users screen - Password Management

4.2.4.2 Password Management

The instrument password can be changed in the **Password Management** area.

The default initial password is:

000000

To change the password:

1. Tap **Users** on the left navigation bar.
2. In the **Old password** field, enter the current password.
 - For a new instrument or a password that has not been changed, enter **000000**.
3. In the **New password** field, enter the new password.
4. Tap **Confirm**.
5. The system saves the new password.

The password must use digital numbers. The note on the right side of the screen indicates that the password should contain **6 to 15 digits**.

Caution: Keep the new password in a safe place. If the password is changed and then lost, password recovery is not available from this function.

4.2.4.3 Local Information

The lower part of the Users screen allows the operator to edit four categories of local information:

- **Company**
- **Center**
- **Laboratory**
- **Technician**

Tap the corresponding tab to enter or select information.

This information is used only for local instrument information storage and viewing. It currently has no additional function.

4.2.4.4 Edit Company Information

1. Tap **Users**.
2. Tap **Company**.
3. Enter or select the company title.

4. Confirm or save the information according to the screen prompt.

The right-side note area displays brief guidance, such as:

“Enter or select Company Title”

4.2.4.5 Edit Center, Laboratory, and Technician

To edit other local information:

1. Tap **Users**.
2. Tap **Center, Laboratory, or Technician**.
3. Enter or select the required information.
4. Confirm or save the information according to the screen prompt.

The entered information is stored locally on the instrument.

4.2.4.6 Recommended Operating Notes

- The initial password is **000000**.
- Change the password only when required by the laboratory or distributor.
- Record the new password and keep it secure.
- Use a 6 to 15 digit numeric password.
- Local Company, Center, Laboratory, and Technician information is for instrument display and viewing only.

4.2.4.7 Troubleshooting

Situation	Possible Cause	Recommended Action
The password cannot be changed.	The old password may be incorrect.	Enter the current password again. If the password has never been changed, enter 000000 as the old password.
The new password is not accepted.	The password format may be incorrect.	Use a numeric password with 6 to 15 digits.
The operator forgot the changed password.	The password was changed and not recorded.	Contact the distributor or service provider. Password recovery is not available from this function.
Local information has no visible effect on	The information is for local storage and	Use this section only for local instrument information management.

Situation	Possible Cause	Recommended Action
analysis.	viewing only.	



Warning: Never clean lenses with alcohol, ordinary tissues, cleaning paper, toilet paper, cotton wool or hand towels, which will scratch the lens surface.



Warning: Do not use volatile materials such as benzene or thinner for cleaning. Use of such materials will discolor and/or damage key device surfaces. Do not use organic solvents. Also, please do not scrape strongly during cleaning.



Warning: You should be extremely careful not to damage the objectives as they are a very sensitive part of the system.



Warning: Please don't let any solvent, lip and erosive material touch the equipment.



Warning: Alternating current power must be well earthed. Alternating current power must be stable. It is prohibited to use power source with high-power equipment together, and it's better to equip with voltage-regulated power supply.

If there is smoke, abnormal smell or strange sound of equipment, cut off the power at once, and contact with the distributor.

When power line is inserted or pulled out, you must hold plug itself, not the power line.



Warning: Avoid cleaning with hard or sharp objects.



Warning: Do no wash the device with water.



Warning: Be sure that the switch is in off ("0") before connecting to the equipment to the power source.



iSperm

www.isperm.com